




THD BANDY

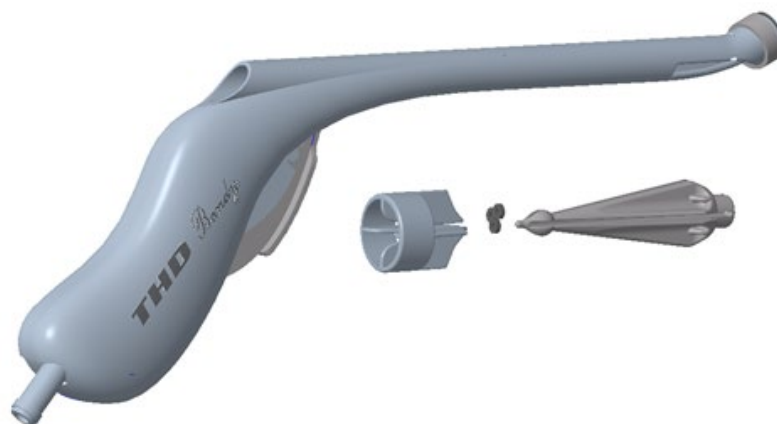

Use

THD Bandy (ref. 800183) is a Class IIa medical device in accordance with the Medical Device Directive 93/42/EEC consisting of a haemorrhoidal binder, a disposable medical device together with a cone and a pusher. THD Bandy is designed for the treatment of degree 2nd hemorrhoidal disease.

General Description

THD Bandy (ref. 800183) is a ready-to-use medical device consisting of an elastic binder together with a cone and a pusher.

The THD Bandy medical device is a disposable suction elastic binder equipped with a system for loading elastic rings supplied separately.

THD Bandy connects to a suction system, to control the quantity of aspirated tissue, around which the elastic band will be released by pulling the trigger.

The special design of the handle allows use with one hand as well as a fast and secure elastic release.

Features

- Made of a medical grade Polycarbonate
- Complies with the biocompatibility requirements of the ISO 10993 standard
- Supplied with disposable rings
- Loading with one hand using loading cone, included
- Ready to use
- Lightweight for easy handling
- CE Class IIa
- Individually packaged
- Latex Free



Advantages

- Can be used with one hand
- Easy loading of the elastic

Dimensional Features THD Bandy (ref. 800183)	
Length except handle (A)	150 mm (5.90 in)
Internal opening diameter (B)	10 mm (0.39 in)
Suction connection diameter (C)	8 mm (0.31 in)

THD Bandy			
REF	Q.ty	Description	Risk Class
800183	25	THD Bandy (medical device) composed of: Disposable suction ligator Cones for loading Ring loader	IIA

Accessoires compatibles	
Product	Description
Rubber band ligator	-Material: Isoprene -Inner diameter: 1.9mm -Outer diameter: 4.5mm