



THD DERM A942 HIGH FREQ.



Use

The Derm 942 High Frequency Desiccator unit is perfect for treatment of anal dysplasia and other High Resolution Anoscopy procedures.

General Description

The Derm 942 High Frequency Desiccator is a low-power electrosurgical unit designed for the office setting and needs low-maintenance level.

The power can be adjusted up to 40 watts and has 1/10th watt increments in the first 10 watts to allow for usage in a wide range of specialties.

Autoclavable Hand Piece with power-switching buttons (800211) and Stainless Steel Stand for mounting on cart are included.

Technical Specifications	
Height	18.8 cm
Width	22.8 cm
Depth	10.5 cm
Weight	< 2.6 kgs
Line Voltage	100-240 VAC
Line Frequency	50-60 Hz
Output Frequency	328-429 kHz



THD SMOKE EVACUATOR SE02-220 ACCESSORIES

THD A902 REUSABLE HANDPIECE

P/N 800211



The Reusable Hand piece is an autoclavable 3-button hand piece intended to be used with the Derm 942 High Frequency Desiccator (800231). Use the blue button for activating the monopolar output. Use the up-down buttons to remotely adjust the power settings of the Derm 942 High Frequency Desiccator (800231). Hand Piece is steam autoclavable.

THD 13-5300 REUSABLE GRINDING

P/N 800209



The Derm 942 High Frequency Desiccator (800231) can be used with or without a grounding cable and pad. To avoid unforeseen shocks to provider or patient, the Reusable Grounding Cable can be plugged into the Patient Plate port of the Derm 942 High Frequency Desiccator and used in conjunction to the Disposable Grounding Pads (800210). This item is tested for 100 uses. Item sold separately.

THD GBS-DB1031A DISP. GROUND.

P/N 800210




The Disposable Grounding Pad is intended to be used with the Reusable Grounding Cable (800209) to ground the patient without physical contact from the provider. This pad works to disperse the energy from the Derm 942 High Frequency Desiccator (800231) and return it safely to the machine without any painful shocks to provider or patient. Item sold Separately.



THD 14-4010HB SMOKE EVAC. FUS. P/N 800212



For simultaneous evacuation of smoke while using the Derm 942 High Frequency Desiccator (800231), attach the Smoke-Evac Fusion High Frequency Desiccator Tubing Kit to the tip of the Autoclavable 3-Button Hand Piece (800211) and with Smoke-Evac Fusion Probes. Attach 9mm connector to Smoke-Evacuator (800181) using the corresponding port.

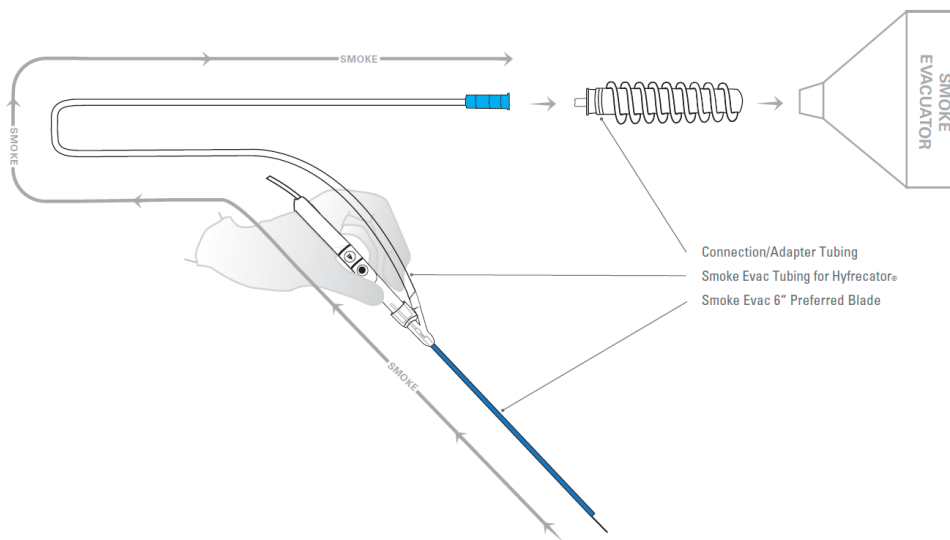

THD SMOKE EVACUATOR SE02-220		
Q.ty	Code	Description
1	800231	THD DERM A942 HIGH FREQ.
1	800211	THD A902 REUSABLE HANDPIECE
1	800209	THD 13-5300 REUSABLE GRINDING
1	800210	THD GBS-DB1031A DISP. GROUND.
1	800212	THD 14-4010HB SMOKE EVAC. FUS.